

Reliable and **complete protection** for safe **disposal** of **cytostatic agents**



the safety system

Emergency cleaning kit

BERNER SpillKit^{Extra} **XP**_{Protection}

Keeping emergencies safely under control

In case of accidents involving carcinogenic substances, it is extremely important to have reliable protective and cleaning equipment at your fingertips. That is why the **BERNER** SpillKit^{Extra} **XP**_{Protection} emergency cleaning kit provides everything you need for safe management and disposal of spilt cytostatic agents. It contains protective clothing, tools, instructions for handling hazardous situations, an accident report form and much more, leaving you perfectly equipped to deal with accidents involving cytostatic agents.


Complete protection certified to the highest standards



You can never take too many precautions for your personal safety. That is why we place the greatest emphasis on obtaining the most important certificates for all of the relevant products. We also have production reviewed by the competent employers' liability insurance association and TÜV on a regular basis. After all, your health is at stake!

All you need for your personal safety

If an accident happens, you won't have to search for long. For maximum safety, we have packed the contents of **BERNER** SpillKit^{Extra} **XP**_{Protection} so that you find them in the order in which you need them. Just below the user information, you will find your

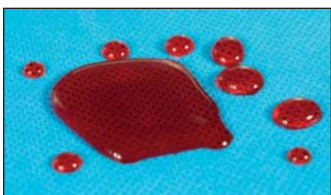

Protective clothing with a coating that is impervious to liquids.

personal protective equipment consisting of a hooded overall, overshoes, 2 pairs of gloves, a respirator mask and safety goggles, all made from safety materials. The reliability of the barrier materials has been proven in comprehensive permeation tests with diverse chemicals, so you will be completely protected from head to toe.

Clear, convenient and easy to use

BERNER SpillKit^{Extra} **XP**_{Protection} comes in two different sizes and has everything you need for dealing with any kind of accident. With its particularly large scoop, sweeper, tongs, distilled water, highly absorbent ChemoSorb pad and tearproof ISYSOFT cloths, you are perfectly equipped to deal with contamination in both liquid and powder form and broken glass. Its clever design and comprehensive documentation ensure that it is used correctly, and the 10-point operating instructions in picture form, which are universally understood, make it particularly user-friendly and safe.



The contents of **BERNER** SpillKit^{Extra} **XP**_{Protection} comply with QuapoS3 [Quality Standards of the Oncology Pharmacy Service].

Order your emergency cleaning kit for cytostatic agents now

Item/size	Item no.
BERNER SpillKit ^{Extra} XP _{Protection} / M	BI-4104
BERNER SpillKit ^{Extra} XP _{Protection} / XL	BI-4114

The following items are included:

Item/size	Quantity/no. of pairs
Protective overall / M or XL	1
Protective gloves / M or XL	1
Protective gloves, yellow / M or XL	1
Protective overshoes / universal	1
FFP3 respirator mask / universal	1
Safety goggles / universal	1
Marker pen	1
Warning sign	1
ChemoSorb pad	1
ISYSOFT cloths	4
Scoop	1
Sweeper	1
Tongs	1
Distilled water	500 ml
Waste bag, blue	1
Special waste bag, white	1
Cable ties	2
User information	1
Accident report	1

Berner International GmbH owns the copyright for the contents of this brochure. Subject to technical and production-related changes.

CE 0299



Made in Germany

GB - 12 0373 A



BERNER INTERNATIONAL GMBH
Mühlenkamp 6 • D-25337 Elmshorn
Postfach 245 • D-25302 Elmshorn
Germany

☎ +49/(0)4121/4356-0
☎ +49/(0)4121/4356-20
✉ info@berner-international.de
🌐 www.berner-international.de



the safety system